

# Solo PIV

## Nd:YAG Laser Systems



### Reliable Operation

- Thermally compensated resonator assures stable operation.
- Rugged I-Beam resonator design on Solo 120XT and Solo 200XT
- Requires minimal maintenance, increasing system up-time.
- Field-proven reliability permits users to concentrate on their applications, rather than on system upkeep.

### Exceptional Performance

- Superior, proven design provides stable, high-energy output with excellent beam quality and pulse-to-pulse stability.
- Compact resonator design provides excellent beam pointing and energy stability.
- Predictable, high performance ensures that your work gets done faster.

**S**olo PIV is a compact, dual laser-head system designed to provide a highly stable green light source for Particle Image Velocimetry (PIV) applications. It is ideally suited for most liquid and many air-based PIV experiments, and its small size provides excellent flexibility in setting-up such experiments.

### Features

- **Small laser head requires minimum space**
- **Single power supply simplifies setup and enhances mobility**
- **High output energy**
  - 15 - 200 mJ at 532 nm
- **Highly flexible design with repetition rates**
  - From 1 to 15, or 30 Hz, depending on model selected
- **Operating convenience provided through multiple triggering capabilities**
  - Continuous internal trigger
  - External TTL trigger
  - Single input pulse activating laser lamp and Q-switch
  - Separate pulses to control lamp & Q-switch independently for precise laser pulse timing control
- **Easy set up:**
  - Single power supply features internal, closed-loop cooling system
  - Operates on 95-240 VAC single phase source
- **Convenient operation made possible with:**
  - Remote positioning of a single power supply - saves valuable lab space
  - Local control panel on power supply with all system controls, including optional optical attenuator



**NEW WAVE**  
RESEARCH™

## Solo PIV Products

		Solo I-15	Solo II-15	Solo II-30	Solo III-15	Solo 120XT	Solo 200XT
<b>Repetition Rate (Hz)</b>		15	15	30	15	15	15
<b>Energy<sup>1,2</sup> (mJ)</b>	532 nm	15	30	30	50	120	200
	266 nm	NA	NA	NA	NA	20	30
<b>Energy Stability<sup>2</sup> (±%)</b>	532 nm	4	4	4	4	4	4
	266 nm	NA	NA	NA	NA	9	9
<b>Beam Diameter (mm)</b>		3	3	3	4	5	6
<b>Pulse Width<sup>3</sup> (ns)</b>		3-5	3-5	3-5	3-5	3-5	3-5
<b>Divergence<sup>4</sup> (mrad)</b>		< 3	<3	<3	<4	<3	<4
<b>Beam Pointing (urad)</b>		<100	<100	<100	<100	<100	<100
<b>Jitter (±ns)<sup>5</sup></b>		1	1	1	1	1	1

- Optical losses due to optional attenuator will reduce maximum energy by 10%.
- Energy and pulse-to-pulse stability for 98% of shots after 30 minute warm up.
- Full width half maximum.
- Full angle for 86% of the energy, at 1/e<sup>2</sup> point.
- From Q-switch synch out to light pulse, 98% of 1000 shots.

## Physical Characteristics


	Laser Head*			Power Supply		
	Solo I, II, III	Solo 120XT	Solo 200XT	Solo I, II, III	Solo 120XT	Solo 200XT
<b>Length</b>	13.775"/350 mm	22.4"/569mm	24.4"/620 mm	18.15" / 461 mm	19.0"/483 mm	21.2" / 538 mm
<b>Width</b>	7.0"/178 mm	9.24"/235 mm	9.24"/235 mm	7.77" / 194 mm	8.6"/218 mm	10.6" / 269 mm
<b>Height</b>	3.187"/81 mm	4.86"/123 mm	4.86"/123 mm	14.32" / 363 mm	15.0"/381 mm	15.16" / 385 mm
<b>Weight</b>	10 lbs./4.5 kg	34 lbs./15.5 kg	40 lbs./18.2 kg	48 lbs. / 22 kg	53 lbs/24 kg	53 lbs. / 24 kg
<b>Length Umbilical</b>	8 ft / 2.4 m	10 ft / 3 m	10 ft / 3 m			

\* Width and height include mounting plate

## Operating Requirements

Temperature	50° - 86° F (10° - 30° C)	
Relative Humidity	20—80% non-condensing	
Voltage	95—240 V, 50/60 Hz	
Power	Solo I, II, III	15 Hz-800 watts; 30 Hz-900 watts
	Solo 120XT	1000 watts
	Solo 200XT	1500 watts

Specifications and product offering subject to change without notice.

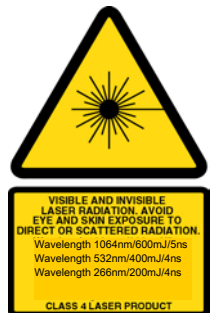

**USA**  
New Wave™ Research, Inc.  
48660 Kato Road  
Fremont CA. 94538-7339  
Tel: 510-249-1550  
Tel: 800-566-1743  
Fax: 510-249-1551  
Email: [Lasers@new-wave.com](mailto:Lasers@new-wave.com)

**Japan**  
New Wave™ Research, KK  
5F Chojiya Building, 1-36-4,  
Shinjuku-ku, Shinjuku  
Tokyo, 160-0022 Japan  
Tel: +81-3-3351-0131  
Fax: +81-3-3351-0121  
Email: [NewWaveKK@new-wave.com](mailto:NewWaveKK@new-wave.com)

**Taiwan**  
New Wave™ Research G. C. Co., Ltd.  
2Fl., No. 58-3, XingShan Rd.,  
Neihu Dist.,  
Taipei, Taiwan 114, R.O.C.  
Tel: 886-2-8792-7585  
Fax: 886-2-8792-7584  
Email: [NewWaveGC@new-wave.com](mailto:NewWaveGC@new-wave.com)

**Europe**  
New Wave™ Research Co., Ltd.  
8 Avro Court, Emine Business Park  
Huntingdon, Cambs  
England PE29 6XS, UK  
Tel: 44-(0)1480 456 566  
Fax: 44-(0)1480 456 545  
Email: [NewWaveEU@new-wave.com](mailto:NewWaveEU@new-wave.com)

**China**  
New Wave Research (China) Co., Ltd.  
Rm. 1102, 11/F, Jiaying Mansion,  
No. 877, Dongfang Rd.  
Shanghai, 200122, China  
Tel: 86-21-5058-7785  
Fax: 86-21-5058-7786  
Email: [NewWaveCN@new-wave.com](mailto:NewWaveCN@new-wave.com)